

MOLEMAXPLUS

DIGITAL DERMOSCOPY SOFTWARE

Enhance patient experience and cost efficiency

MoleMax Patient Report Report Date: XX XX XXXX

Checklist Name: _____ Clinic Name: _____

Patient Details: Name: XXXXXXXXXXXXXXXXXX Age: XX Sex: M Patient Record #: XXXXXXXXXXXXXXXX Last Visit Date: XX XX XXXX

Lesion Clinical Observation: Location: Right lower leg Date: XX XX XXXX

Prescribed Diagnosis: _____ Requested Histology: _____

Suggested Management: _____ Exclusion: AA AA AAAA

Histopathology Report (if available): _____



Final Diagnosis: Melanoma

Comments: _____

Next Appointment Date: XX XX XXXX

Clinic: _____ Phone: _____
 Clinic Address: _____ Email: _____
 Clinic Phone Number: _____

Page 1 of 1

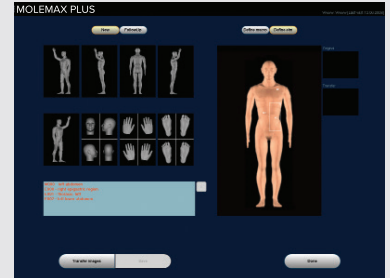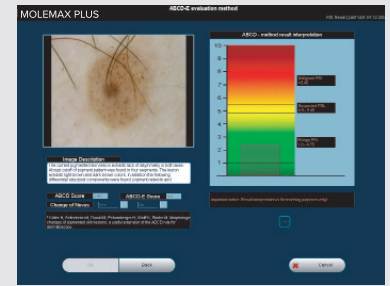

MoleMax Plus is a software package for dermoscopy which was specially developed for diagnostic processing with images captured from any dermatoscope using a compatible smart phone or digital cameras. Simple clinical and overview imaging makes trending and monitoring an easy task. MoleMax Plus runs on all the latest Windows operating systems.

- ▶ FREE 30 Day Trial
- ▶ Link to MoleMax GO app for easy record entry on the go
- ▶ Localise patient images easily using the Body Map
- ▶ Monitor the development of a lesion over time
- ▶ An All Skin diagnostic atlas with clinical references
- ▶ Integrated diagnostic aids such as the widely used ABCD rule & the Seven Point rule
- ▶ Compatible with a wide range of dermatoscopes with smartphone connectivity
- ▶ Free to use with MoleMax Plus or use as a standalone setup
- ▶ Links with your Practice Management Software (PMS)

MOLEMAXGO

DIGITAL DERMOSCOPY PHONE APP

Automatically link images taken on your smartphone or tablet with MoleMax Plus and your patient management software. Available for iOS and Android devices.

- ▶ Compatible with Third Party Dermatoscopes
- ▶ Easy to Use Patient Report Function
- ▶ Free Local Secure Storage
- ▶ Choose from a one off fee or a basic subscription



MOLEMAXPLUS

DIGITAL DERMOSCOPY SOFTWARE

Which MoleMax System is best for you?

MoleMax Range Comparison Matrix*

Software Features	MoleMax HD PRO	MoleMax HD	MoleMax HD Mini/Modular	SkinDOC	MoleMax Plus	PhotoMAX
All Skin Diagnosis Reference Library	✓	✓	✓	✓	✓	✓
Practice Management System Reporting Output	✓	✓	✓	✓	✓	
Comparative Screening / Follow Up	✓	✓	✓	✓	✓	✓
One Click Automatic BackUp	✓	✓	✓	✓	✓	✓
Monitoring / Trending Follow Up	✓	✓	✓	✓	✓	✓
MoleMax Go Integrated Phone App(Android / IOS Apple)	✓	✓	✓	✓	✓	
Mole Scoring Expertizer Plus Scoring Aid	✓	✓	✓	✓	✓	
Lesion Measurement Aid	✓	✓	✓	✓	✓	
Total Body Mapping (Import Images)	✓	✓	✓	✓	✓	✓
Total Body Mapping (Live Images)	✓	✓	✓		✓ Via Mobile	
Live Video High Definition (HD) Camera Imaging	✓	✓	✓	✓	✓ Via Mobile	
Live Digital SLR Camera Imaging	✓	✓	✓			
Thumbnail Follow Up	✓	✓	✓			
Clinical Audit and Quality Control	✓	✓	✓			
CosMax Cosmetic Simulation Module	✓	✓	✓			
Education and Research Module	✓	✓	✓			
Software Operation via Mobile Tablet	✓	✓	✓			
Touch Screen Operation	✓	✓	✓	✓	✓	✓
SQL Database Integration	✓	✓	✓	✓	✓	✓
Trichoscan Module	✓	✓	✓	✓	✓	
Additional User Licence on networked laptop	✓	✓	✓	✓	✓	✓
Multi-User Licence	✓	✓	✓	✓	✓	✓

Hardware Features	MoleMax HD PRO	MoleMax HD	MoleMax HD Mini/Modular	SkinDOC	MoleMax Plus	PhotoMAX
Camera Docking Station	✓	✓	✓	✓		
HD Camera - Micro Imaging (30x Only)				✓		
HD Camera - Micro & Macro Imaging (Variable Magnification)	✓	✓	✓			
High Resolution Digital SLR Camera for Macro Imaging	✓	✓	✓			
Polarised Adapter (20/30/40x) for Dermoscopic Imaging	✓	✓	✓			
Polarised Adapter (60/80/100x) for Dermoscopic Imaging	✓	✓	✓			
Non-Polarised Adapter (20/30/40x) for Dermoscopic Imaging	✓	✓	✓			
Non-Polarised Adapter (60/80/100x) for Dermoscopic Imaging	✓	✓	✓			
Networked Windows Surface PRO Tablet or Laptop	✓	✓	✓			
USB 3.0 Back Up Hardware	✓	✓	✓	✓ Mini PC Opt	✓ Mini PC Opt	✓ Mini PC Opt
2nd Back Up Storage Drive	✓	✓	✓			
Full Spec Intel i7 PC with High Speed Hard Drive	✓	✓	✓	✓ Mini PC Opt	✓ Mini PC Opt	✓ Mini PC Opt
Full Spec Intel i5 PC with High Speed Hard Drive		✓	✓	✓ Mini PC Opt	✓ Mini PC Opt	✓ Mini PC Opt
Monitor	✓ 2 x Touch 23 Inch	✓ 2 x 23 Inch	✓ 1 or 2 x 23 Inch	✓	✓	✓
Keyboard	✓	✓	✓	✓	✓	✓
Mouse	✓	✓	✓	✓	✓	✓
Trolley Design for Mobility with Monitor Mounting	✓	✓	✓			
Manual Total Body Mapping Stand	✓	✓	✓			
Automatic Total Body Mapping Stand	✓ Coming Soon	✓ Coming Soon	✓ Coming Soon	✓ Coming Soon		

FEATURE KEY:	Included	✓
	Optional	✓
	Not Available	

* As at 1/6/2023

MOLEMAXPLUS

DIGITAL DERMOSCOPY SOFTWARE

MOLEMAX PLUS Software Package

Free 30 Day Trial

MoleMax Plus is an entry level software package for digital dermoscopy which was specially developed for skin screening with digital cameras and smartphones. Simple clinical and overview imaging makes trending and monitoring an easy task.


MoleMax Plus is a medical program package which features:

- ▶ Localise patient images easily using the Body Map
- ▶ Monitor the development of a lesion over time
- ▶ The All Skin diagnostic atlas with clinical references
- ▶ Integrated diagnostic aids such as the widely used ABCD rule & the Seven Point rule
- ▶ Links with your Practice Management Software (PMS)


All Skin Diagnosis Library: Diagnosis can be selected by names and images showing various skin conditions. Customised Patient Reports: Print images, patient data and history along with results from the scoring module. Backup and Export: Easy backup of patient history, transfer or share databases between PCs in your clinic.

Minimum System Requirements:


Windows 7 or above, Intel i3 or above,
Graphic Display Adapter Capable of
1280x1024 True Colour (24 or 32 bit), 100 MB
free Hard Disk Drive space (800 MB with
Mole Atlas), CD-ROM Drive.

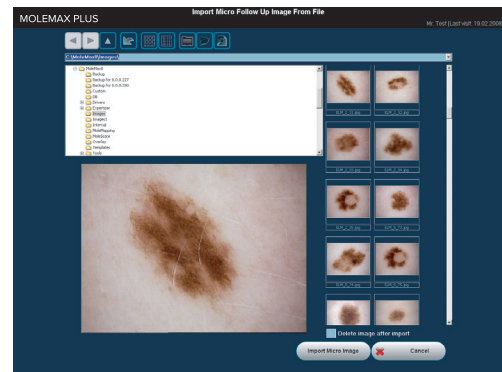


Import Images Taken from Multiple Digital Sources

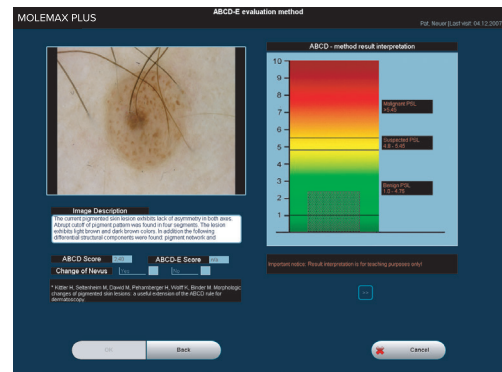


ABCD and 7 Point Rule Scoring and Measurement

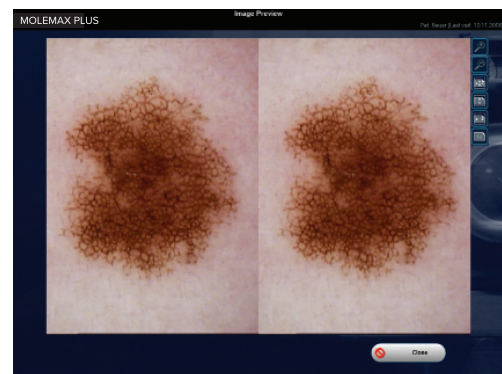

Comparative Imaging and Follow-up

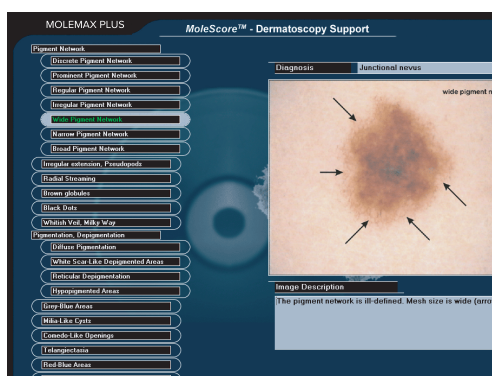

Extensive Image Library of All Skin Diagnoses


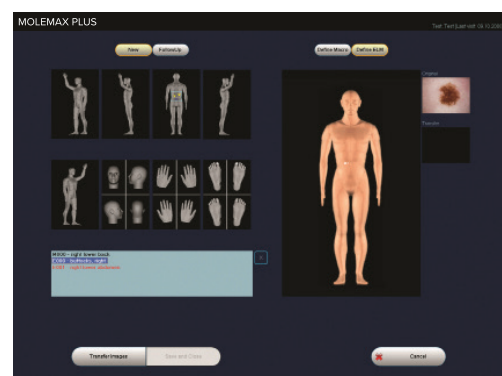


Body Map for easy lesion localisation

MOLEMAXPLUS

DIGITAL DERMOSCOPY SOFTWARE

Software Features

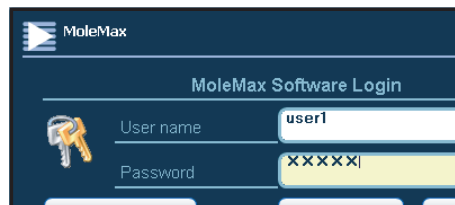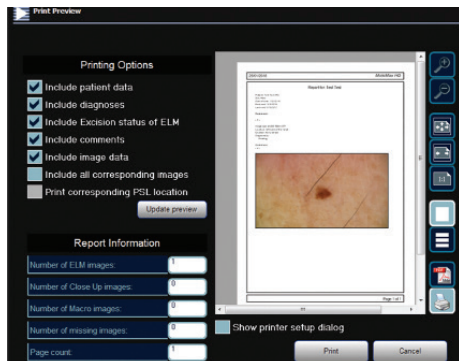

Multi-User Environment:

Our software is specially designed to operate in a multiuser environment to maintain audit paths and patient data security. The Administrator may tailor access rights for each user, to limit read, write and delete access.


Data Back Up & Export:

The MoleMax and SkinDOC systems runs perfectly within most clinic networks. Data can be stored and backed up on a clinics' server. Where a spare server is not available, the MoleMax is equipped with ports for external backup. Patient images and files may also be exported from the MoleMax system into practice management software.




Real Time/Overlay Follow Up:

A quick comparison in large scale is one of the most important features Doctors look for in Skin Imaging software. A follow-up image may be taken with the previous image in view on one side of the screen. Alternatively, an outline may be generated of the original image to allow for easier follow-up. Images over time may also be compared side by side or over each other to easily assist in recognising lesion alterations.

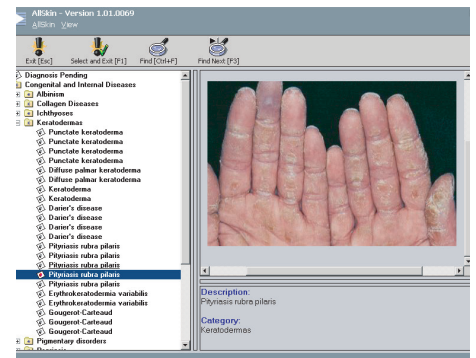


All Skin Library:

The All Skin module is a diagnosis library filled with images and descriptions which give the user the ability to browse and search for certain skin conditions. The skin conditions vary from melanoma to actinic keratosis, basal cell carcinomas and many more. The skin library features hundreds of them, including those not related to skin cancer.

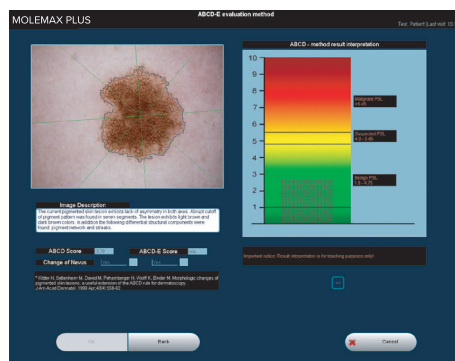

Expertizer Plus Scoring Aid:

This module contains a histopathology database of diagnosed pigmented lesions taken on the MoleMax system. These diagnosed images can be used for comparison and scoring purposes. This module also includes image analysis functions which will provide the user with automated calculations for various parameters such as diameter, circumference, and area measurements. The module can also auto-measure and compare any two lesions taken over time to analyse pattern and size changes. The aid-to-diagnose module has in-built scoring tools which guide the user to assess risk of a lesion with a score.

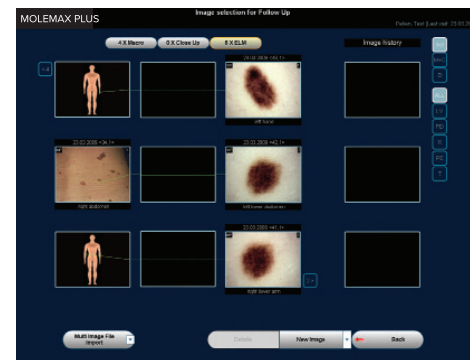

Trending and Monitoring:

Macro Trending allows for Micro and Histopathology images to be taken over time with comments. A complete patient history can be kept in chronological order every image taken and every comment added to the patient file.